



**EunSung**

Innovative Technology For your Medical Aesthetics Needs





# EunSung



# 3-MAX<sup>plus</sup>™

**The Optimum Solution for Cellulite and Body Shaping Treatment.**



## Ultrasonic cavitation, RF with Vacuum and expert RF system

Clinical power shaper combines Cavitation, RF and Vacuum technologies. 3 RF applicators are updated in -3Max plus system. The applicators combine Multi-polar radiofrequency and LED technology.

### Specification

Radiofrequency	Multi-polar (1MHz), Bi-polar (2MHz)
Ultrasonic frequency	36 KHz
Vacuum Range	370 mmHg
Input Power	AC 110~220V, 60/50HZ
Weight	53Kg
Dimension ( WxDxH)	406 x 505 x1,252 (mm)

### • Type of handpiece



#### RF body applicator

- Specialized dynamic vibration function
- Easy to use technical motions
- Updated multi-polar RF
- LED technology (● Red)




#### Cavitation body applicator

- Powerful cavitation effect from dual generators
- Effective treatment by 4 types of pulse mode
- Guaranteed treatment results
- Operating signal light function



#### RF facial applicator

- Dramatic result for facial lifting
- Dynamic & Strong vibration function
- Effective facial lymphatic drainage
- 4 types of vacuum pulse

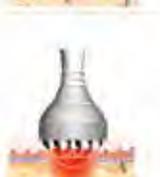

#### RF body & Facial Expect Applicators

- Stable multi-polar RF system (1MHz, 2MHz)
- Specialized body, arm & leg and facial applicators
- LED technology (Large & Medium: ● Red / Small: ● Blue)



Large (Body)



Medium (Arm&Leg)



Small (Facial)



# POWER SHAPE 2™

## Multi-functional body shaping system

Powershape2 is an intergraded system for body contouring and face lifting combined with Multi & Bi-polar RF, Thermal RF, Vacuum and Low Level Laser.

### ▪ Structure of handpiece


#### Roll Handpiece

Bi-polar Vacuum-assisted motorized roll handpiece enables both laser & RF energy precisely target the areas as to where it can have maximum effect.



#### Multi-polar handpiece

Multi-polar RF technology painlessly delivers focused and optimal RF energy to the skin surface and to subcutaneous fat layers.

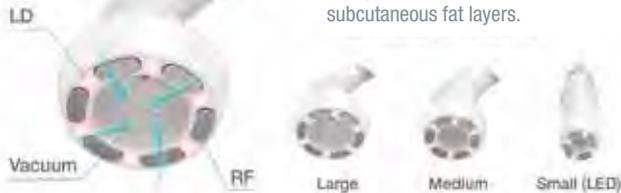

#### Thermal handpiece

Button

Increased thermal energy conductivity that reduces treatment time thanks to the ceramic coated bi-polar tips.



### Specification

Type	Radiofrequency, Vacuum, Low Level Laser
Applications	Roll handpiece, Multi-polar handpiece, Thermal handpiece
Panel	10.2 inch, Digital touch screen
Input Power	AC 100~240V, 60/50Hz
Weight	60Kg
Dimension ( WxDxH)	414 x 523 x1,230 (mm)









**CRISTAL TECHNOLOGY:**

**The widest range of applicators on the market**  
 CRISTAL has 8 different applicators that can treat all parts of the body in men and women. Every year, Deleo's R&D department patents new technologies to offer you more innovative and powerful devices.

**Safety first** CRISTAL is the safest body contouring device on the market with 15 control probes. CRISTAL has been developed in accordance with European medical standards and has been granted the Medical CE mark.

**Protocols yielding better results** CRISTAL offers many pre-recorded treatment protocols that have been studied and developed by our committee of scientific experts.

**No consumables, more profitability** CRISTAL ensures you a high profitability for a fast development of your activity.

**Made in France** CRISTAL is produced in France in Saint-Raphaël, in the Var region. The DELEO know-how and expertise guarantee a high technology, reliable and secure device.



**BEFORE / AFTER RESULTS:**

Before



After



Before



After



**A BROAD CHOICE OF APPLICATORS TO TREAT ALL BODY PARTS:**



A1; small size for chin



Emerald small fat clusters for arms, fold of the back&knees



A2; small fat deposits for knees & chest area



Shapphire all type of fats for abdomen & thighs



A3; medium fatty clusters for multi-body areas



Ruby fatty clusters of significant size for abdomen



Diamond cellulite reduction, all body areas



Amethyst flat areas for horse panties,calf & thighs



TRADEUM Medical Solutions Co.W.L.L. شركة تراديوم للحلول الطبية ذ.م.م

A Part Of Jothen

إحدى شركات الجوثن

www.tradeum.com

1800166

tradeumbyaljothen

tradeum





**CRISTAL Cryolipolysis is a medical device intended to treat localised fat rolls for patients suffering from weight problems and undergoing certain diets. It allows to improve the skin texture and to significantly diminish localized fat rolls.**

### **WHAT IS MEDICAL CRYOLIPOLYSIS ?**

Medical cryolipolysis is a non-invasive technique of adipocyte lysis. The principle of cryolipolysis is to intensively cool a localized fat mass with Peltier modules that go down to  $-9^{\circ}\text{C}$ . High pressure is also applied to the adipocytes by aspiration to induce ischemia, in order to modify their configuration. This change in cellular state from a normal body temperature ( $37^{\circ}\text{C}$ ) to very low or even negative temperatures sends a chemical message of apoptosis to the adipocytes. Results are visible 6 weeks later, but the body will continue to eliminate crystallized adipocytes up to 3 months after the session.

### **Powerful, stable cooling system**

The CRISTAL® cryolipolysis system uses an R134a gas cooler, allowing it to quickly descend to very low temperatures while maintaining a stable and steady cold temperature throughout treatment. This technological innovation guarantees satisfactory results for your patients.



TRADEUM Medical Solutions Co.W.L.L شركة تراميوم للحلول الطبية ذ.م.م

A Part Of  
jothen

إحدى شركات  
الجوثن

[www.tradeum.com](http://www.tradeum.com)

 1800166

 [tradeumbyaljothen](https://www.facebook.com/tradeumbyaljothen)

 [tradeum](https://www.instagram.com/tradeum)

